


FLUSSO™

By Pass Adapter

To maintain a closed circuit



Integrated Swing Valve Technology™

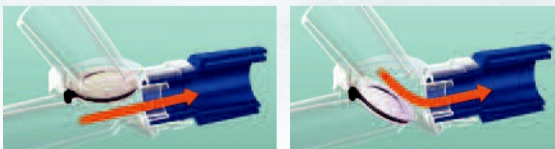

Made in Canada

Flusso™ and Swing Valve Technology™ are trademarks used under licence by
McArthur Medical Sales Inc.

20  18

SILVER WINNER

**MEDICAL DESIGN
EXCELLENCE AWARDS**


US and International Patents Pending



Flusso™ Highlights:

- Proprietary valve with Swing Valve Technology™
- Clear housing with colored valve permits visualization of function of Swing Valve Technology™
 - Integrated swivel at patient connection to reduce torque on endotracheal tube
 - Tethered dust cap keeps ports clean and protects swing valve
- Minimal dead space (7cc)
- Standard ISO 15/22mm connections

Flusso is a by pass adapter that can be used with passive or active humidity systems to maintain a closed circuit. Flusso facilitates the disconnection of the patient from the mechanical ventilator during a circuit change, HME change, or disconnection of patient for transport.

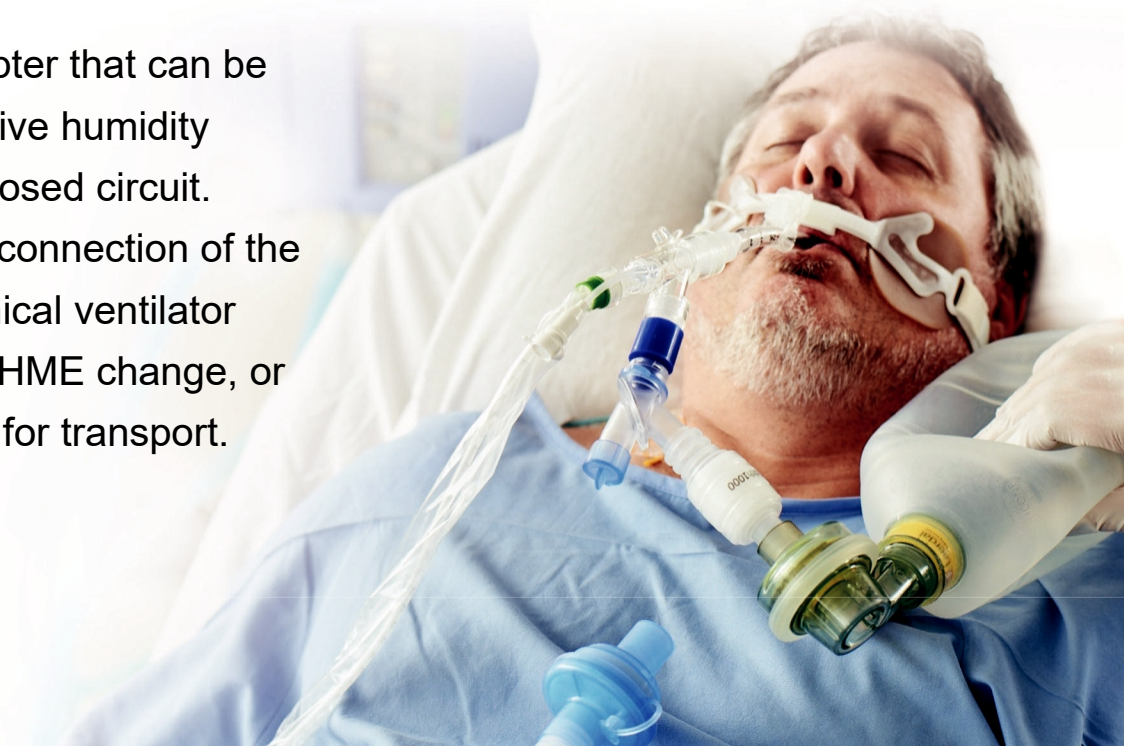

FLUSSO[™]
BY PASS ADAPTER

www.flussobypass.com

To Order: Product Code 420-1978 25/case

Distributed by:

**SOVEREIGN
MEDICAL, INC.**

620H Valley Forge Rd
Hillsborough, NC 27278

919-644-1113

sovmed.com

Manufactured by:
Southmedic Inc.
50 Alliance Blvd.,
Barrie, ON, Canada,
L4M 5K3

20  18
SILVER WINNER
MEDICAL DESIGN
EXCELLENCE AWARDS

10th INNOVATION
PRIZE PARIS-2018
An event by Universal BioTech